This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and the receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.
- Increase the separation between the equipment and the receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

Caution: Any changes or modifications to this device not expressly approved by the manufacturer may void your authority to operate this equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The device has been evaluated to meet general RF exposure requirements. The device can be used under portable exposure conditions without restriction.

£ £ RoHS

M

Product model may change due to improvements without notice.

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RF EMS BEAUTY INSTRUMENT

Lead in Nutrition
Clamping and lifting

Nourishing & Anti-wrinkle
Eliminate acne and whitening

Suitable for eyes, face, neck and body

I INTRODUCTION

Over time, the skin may become wrinkled, discolored, and saggy. This is due to normal sun exposure, natural aging, acne, genetics, etc. Skin care products, such as serums, lotions, and moisturizers, work at the surface level of the skin but generally cannot penetrate deep into the skin and affect collagen at the deeper layers of the skin. the skin.

Radio Frequency (rf) and Electrical Muscle Stimulation (EMS) beauty instruments work with unique, deeper stimulating effects that can stimulate collagen regeneration, firm skin, and improve skin tone and overall appearance skin. Moisturizers and lotions are encouraged to be used with this device as it can increase their effectiveness and ensure smoother operation of the device

This device is recommended for the following conditions:

- softens fine lines and wrinkles
 - for UV-damaged skin to prevent skin aging
 - for acne-prone skin
 - for skin with imperfections against pigment spots
-

For many people, this device can eliminate the need for costly and painful injections aimed at delivering skin products to the deepest levels of the skin. This device works on five fundamental concepts of healthy and painless skin care. These concepts are electrical muscle stimulation (EMS), mesoporation, electroporation, radio frequency (rf) and light emitting diode (LED) light energy.

1. Ems (Electrical Muscle Stimulation)

Ems is used to activate skin cells, stimulate collagen rejuvenation, improve skin elasticity, soften wrinkles and leave skin feeling smooth and healthy.

2. Mesoporation

The origin of mesoporation comes from the combination of mesotherapy and electroporation. Mesotherapy is a transdermal skin delivery system in which micro-needles puncture the skin to allow delivery of active ingredients to deeper layers of the skin

Electroporation is a non-invasive electrical treatment that also transports active ingredients into the dermal layers. Unlike mesotherapy, mesoporation does not require needles.

Mesoporation works by disrupting the surface of the skin, causing the disrupted cells to temporarily open up to create tiny moisture channels or mesopores, which allow the transport of high molecular weight agents into the skin.

3. Electroporation

Electroporation provides a way to transport collagen and hyaluronic acid deeper into the skin. Electroporation provides the means to achieve better results with your skincare products without the need for needles.

4. Radio frequency (RF)


Radio frequency energy uses the natural resistance of tissues in the different layers of the skin to transform RF energy applied to the skin into thermal energy. RF energies can be used for patients of all skin types and allow different depths of penetration allowing for increased collagen contraction and the production of new collagen.

5. Light emitting diode light energy (LED)

The LED light is used to provide painless skin care that plumps the skin by stimulating collagen production and treats mild to moderate acne. The LED light works by sending energy-producing packets of light into the deeper layers of the skin. LED treatments can minimize fine lines and wrinkles, treat sun damage and stretch marks, and reduce redness after more aggressive IPL or more aggressive treatments.

II NAME OF PARTS

- | | |
|-----------------------|---------------------------|
| (1) Charging stand | (2) Charging socket |
| (3) USB charging line | (4) Conductive metal head |
| (5) LED panel | (6) LCD screen |
| (7) Power switch | (8) LED selection key |
| (9) Charge indicator | |

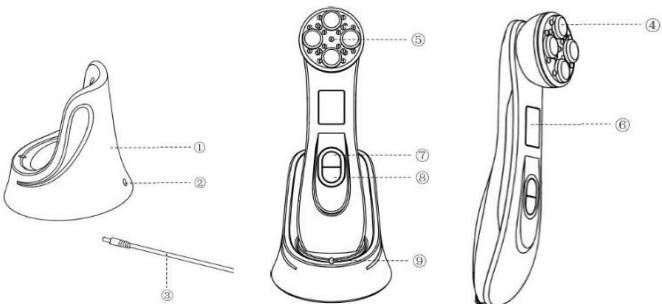III INSTRUCTIONS FOR USEI

1. Charging - Charge the device for 3 hours before first use. To charge, insert the connector of the USB cable into the device and the other end of the cable into the charging socket and connect to the power source.

2. After cleansing the skin, apply your favorite serum to the skin.

3. Press the power switch on the device. The red LED will light up and the device will be at intensity level 1.


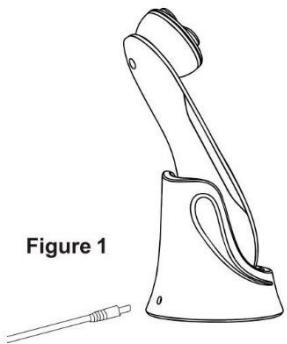
Figure 1

4. Adjust intensity level based on skin type and personal preference

5. Adjust the intensity by pressing the power button, press once for intensity level 1

- Press twice for intensity level 2
- Triple tap for intensity level 3
- Press four times for intensity level 4
- Tap five times for intensity level 5
- Press six times to turn off the device

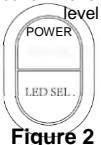6. Slowly move the device up and down the face and from inside to outside. Skin may feel warm during use and this is normal. The device will automatically turn off after 10 minutes of use

It is recommended to use the device 3 times a week for 10 minutes each time. It is also recommended to use a skin serum to keep the skin moist during each treatment with this device.

IV. HOW LED LIGHT PROCESSING WORKS

Mode 1- Red LED Light Treatment Red LED light is used to treat wrinkled skin and improve collagen regeneration and skin elasticity. Press the LED select button once to enter mode 1.

Mode 2- Green LED Light Treatment Green LED light is used to treat dark spots or discolored skin by preventing overproduction of melanin. Press the LED select button twice to enter mode 2.

Mode 3- Blue LED Light Treatment Blue LED light is used to treat acne and oily skin. Press the LED select button three times to enter mode 3.

Mode 4 - Yellow LED light treatment

A yellow LED light is used to lighten the skin. It is suitable for skin with dark spots and dull edges. Press the LED select button 4 times to enter mode 4.

5-Pink Mode LED Light Processing

Pink LED light is used for skin whitening. Press the LED select button 5 times to enter mode 5.

Mode 6 flashing pink LED light treatment

The pink flashing LED light is used to facilitate the permeation of skin base products. Press the LED select button 6 times to enter mode 6. Turn off the LED light

Press the LED select button 7 times to turn off the LED light function.

Correct method of operation Metal head in full contact with the skin

V PARAM

Incorrect method of operation Metallic partial contact skin

N

Name: RF & EMS Beauty Instrument Model: 9902

Input: DC5V 500mA

Battery Capacity: 750mAh

Power: Max 10W

Charging time: 3.5 hoursStandby time: 40 days

Use time: 10 minutes / hour

Unit weight: 220g

Unit size: 165*45*46mm

Gift box: 215*152*85mm

VI SAFETY INSTRUCTIONS

1. Read all safety instructions before using and use it correctly.
2. Follow the user instructions to use this device.
3. Do not lower or redo the device, please consult the distributor for repair.
4. Use the device only for the purpose for which it was designed, as described in this manual.
5. Do not store and use the device in a place with high temperature and humidity such as under the sun or near a tire or in the bathroom, etc..
6. Do not use or place the device in wet conditions, short circuit, electric shock or tire damage.
7. Keep away from children and physically handicapped people Allow use, close supervision is necessary when this device is used by user under 18 years of age.
8. Keep out of reach of children.
9. Pregnant women and people with hypertension are not allowed to use heart disease
10. The person who is excessively sensitive to skin disease is not allowed to use.

WARNING

To avoid malfunction or fire, please stop using in case of unusual heat.

2. Please do not drop or hit the device.
3. Please apply serum constantly during use.
4. Please do not use it if the power cord is damaged.
5. Below conditions are not allowed to use:

Heart disease, fever, communicable disease, cáncer, open cuts, acute illness, hypotension, hypertension, installation of a pacemaker, artificial heart or lung, during medical treatment, skin disease, sunburn, susceptible people, anapetia

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