## User Manual

# MULTIFUNCTIONAL FACIAL SKINCARE MASSAGER

### Cautions when using:

1. Press and hold the power button for 3 seconds to turn it on. It cannot be turned on when charging.
2. When using, hold the sensor strips on both sides of the wand, with the massage head in contact with the skin.
3. Please use it with a cotton pad when exporting the pores and remove the makeup remover or cleansing lotion with a cotton pad. Do not use it with facial cleanser/cleansing oil
4. This product uses micro-current ion technology. If you feel a small tingling sensation this is normal.
5. When using, do not stay in the same place for more than 3 seconds to avoid to skin irritation.

## Functions

## Cleansing Rejuvenation Device

A facial beauty device with galvanic ion, heating, micro-vibration, and red light. This device aids in deep cleansing, rejuvenation, and facial lifting.


\*Galvanic ionic: Using the principle of a mild impulse current between the anode and cathode, it can soften and clean dirt and oil sebum deposits that lodge in the pores. This deep treatment acts as a natural cleanser that helps reduce pore size, and enriched charged oxygen ions activate nutrients to repair damaged skin.

\*Micro Vibration: Using physical vibration, and electric stimulation are both ways to increase elasticity and enhance cellular metabolism.

\*Warm Massage: Promotes absorption of skin-care products, helps makes skin texture tight and keeps the skin more elastic.


\*Red Light: 640nm-650nm red light deep into the dermis, helps stimulate new collagen and elastin growth and also reduces wrinkles and improves skin elasticity.

## Functions
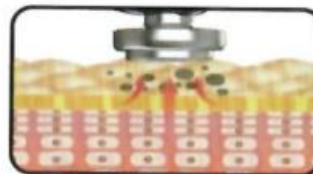
## WARM CLEAN MODE

Warmth helps open the skin's pores, deep cleaning the skin to help electively refine the pores and aid in removing dark skin spots.
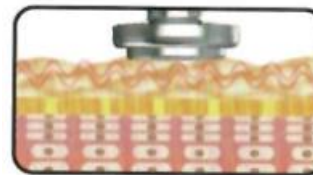
## SKIN LIFTING MODE

Promotes blood circulation, enhances skin elasticity and firming, and creates the best skin condition possible.

## PHOTOTHERAPY MODE

Drive skin-care products deep into the skin, activating the cell's elasticity. Using the red light helps aid in making the skin whiter and smooths out the wrinkles.

## EYE CARE THERAPY MODE

Introduce eye cream or your favorite eye serum adequately, promote eye cream or eye serum absorption, achieve the maximum result of using your favorite eye

## Components

## Use of the Four modes

## Product Parameters:

Product Name:	Photon Rejuvenation Wand
Model No:	W- 618
Product Size:	4.3*16.8cm
Product Weight:	About 90g
Package Contains:	1'device; 1'instruction manual, 1*USB cable
Input:	DC 3.7V/0.6A
Rated Voltage/ Rated Capacity:	3. 7V / 350mAH
Charging Time:	2-3 hours
Intensity:	3 levels
Country Of Origin:	China
Carried Standard:	GB4706.1-2005


Note: In view of the consistent policy of continuous product improvement, the company reserves the right to improve product specifications, the actual product may vary.

## Cautions:

- 1 Clean and rejuvenate, clean the dirt first and then use the cleanser (No need to use the device when cleaning face). Deep cleaning will pull the dirt from deep in the pores to the skin surface, and then wash the dirt with facial cleanser.
- 2.Each mode is timed for 4 minutes. After 4 minutes, it will enter standby mode. After 1 more minute, it will automatically shut down after no action.
- 3.When using, you need to hold the ion sensor strip on both sides of the fuselage and use the massage head to contact the skin to form a circuit. Please keep the skin moist when using otherwise, the device will not work properly.
- 4.If you use it on someone else's face, you need to contact the user's skin to form a circuit, otherwise, the device will not work properly.
- 5.Please use the device according to the length of time specified in the manual. Do not overuse it


**1. Guide to Warm Cleanse Mode** (Requires makeup remover/clean lotion and cotton pad)

This mode can effectively clean the dirt in the pores and make the skin shiny  
Frequency of use: 2-3 times a week. Sensitive skin users should appropriately reduce the frequency of use. Use time: no more than 4 minutes/time, and 1-2 minutes/time for sensitive skin.

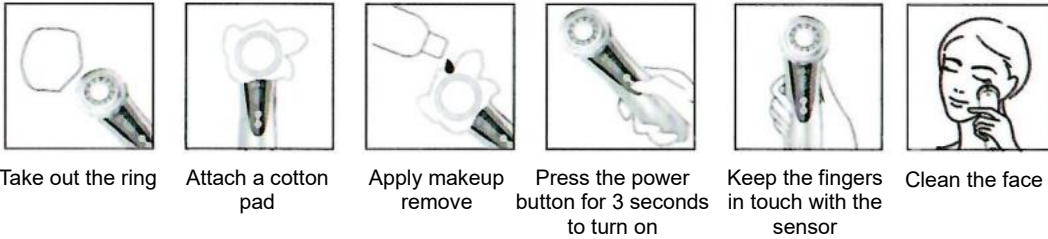
**2. Guide to Phototherapy mode** (requires toner, lotion, essence)

This mode enhances cell activity, diminishes wrinkles and helps skin-care products deeply moisturize the skin. Frequency of use: Can be used once a day in the morning and evening, sensitive skin users should appropriately reduce the frequency of use.  
Use time: no more than 4 minutes/time, 1 -2 minutes / time for sensitive skin, please stop using the device when skin-care products be absorbed by skin.

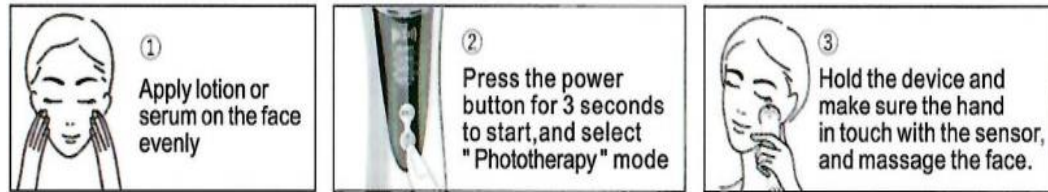**3. Skin Lifting mode** (use with a mask)

This mode gently revives and relaxes the skin, helping to restore the skin's luster and vitality.  
Frequency of use daily

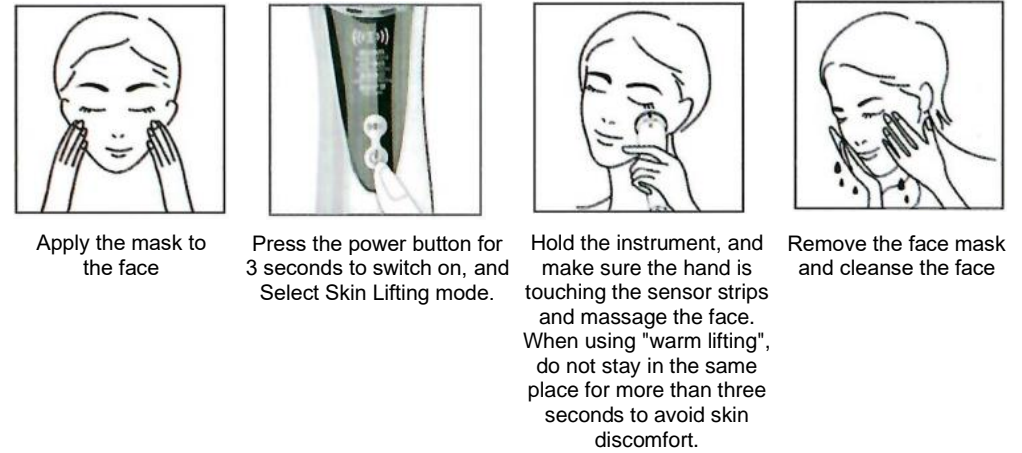**4. Guide to the eye care mode** (use with eye cream)

This mode can fully introduce eye cream or eye essence, promote the absorption of eye cream or eye essence.  
Frequency of use: Can be used once a day in the morning and evening, sensitive skin users should appropriately reduce the frequency of use.

Use time: no more than 4 minutes/time, 1-2 minutes / time for sensitive skin. Do not apply more eye cream or eye essence after absorption.

①

Eye care products (such as eye cream, eye mask, etc.) are evenly applied around the eyes.

②

Press and hold the power button for 3 seconds to switch on, and select "Eye Care" mode.

③

Hold the instrument, touch the sensor strip on both sides of the instrument with your hands, use a point-like pressing technique and massage the head gently from the inner corner of your eyes to the end of your eyes for point-pressing massage. Red light for eye care is off by default and can be turned on automatically by pressing the button as required.

Recommended Practices

Chin

Place the device horizontally at the lower center of the chin, with curved motion and in rhythm for 1-2 seconds, glide the device from left to right and reverse.

Face contouring

Place the device on the cheek. Moving towards the cheekbone, glide the device from bottom to top in rhythm for 1-2 seconds.

Cheeks  
There are two parts to the cheeks [1] & [2] Glide the device from left to right and reverse. Reverse the process alternately

Place the device horizontally at the corner of the lips. Moving towards the cheekbone, glide the device from bottom to top in a rhythm of 1,2 seconds,

Place the device vertically to the corner of the lips. In a rhythm of 1,2 seconds, glide the device horizontally from the corner of the lips toward the ear.

Eye

With eye cream or eye essence, use the point-like pressing technique to gently massage the head from the inner corner of the eye to the end of the eye. Stay at each point for 2 seconds to relieve fatigue and promote absorption.

Recommended Uses Maintenance Aspects

**Morning refreshing**  
Awaken skin + soothing edema +  
smoothen makeup

- Clean face
- Apply toner, essence and lotion
- Phototherapy skin rejuvenation mode
- Makeup

**Evening repair**  
Promoting Skin Absorption +  
Strong Moisture and Nutrition + Overnight Repair

- Deep Cleaning (2-3 times a week in skin clean and rejuvenation mode).
- Clean face
- Apply toner, essence and lotion (phototherapy skin rejuvenation mode)
- Masking introduction (Tempering lifting mode

Please do not wash the wand body  
The instrument will malfunction, please use a dry cloth or small soft cloth.

Cleaning protective cover  
If the cover gets dirty wipe with a clean soft cloth, rinse with water.

Do not wipe the main wand body with a corrosive solution Please do not use alcohol, water, or make-up remover wipe, This may cause product failure or parts cracking, discoloration, etc.

Safety Precautions

1. Not suitable for facial plastic surgery and subcutaneous implants; Please use under the guidance of your doctor if pregnant.
2. Do not wash the device. You should not contact the Photon head with metal objects; To clean use a cleansing wet wipe.
3. Do not put it in the bathroom or wet place, this could cause a short circuit.
4. This device is suitable for the face and neck, however, please avoid pustules and inflammatory acne.

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